ASAHI Conquest Pro /

Confianza Pro

PCI Guide Wire

SYMBOLS	1
English [EN]	4
Dansk [DA]	7
Nederlands [NL]	10
Suomi [FI]	13
Français [FR]	16
Deutsch [DE]	19
Ελληνικά [EL]	22
Magyar [HU]	25
Italiano [IT]	28
Polski [PL]	31
Português [PT]	34
Español [ES]	37
Svenska [SV]	40
Česky [CS]	43
Norsk [NO]	46
Slovenský [SK]	49
Eesti keel [ET]	52
Latviešu [LV]	55
Lietuvių kalba [LT]	58
Türkçe [TR]	61
Русский [RU]	64
Български [BG]	67
Hrvatski [HR]	70
Română [RO]	73
Slovenščina [SL]	76
Malti [MT]	79

SYMBOLS

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PCI Guide Wire

INSTRUCTIONS FOR USE

Read these instructions carefully before using this guide wire and observe the Indications for Use, Contraindications, Warnings, Precautions, and How to Use sections in this Instructions for Use.

Failure to do so may result in complications, including serious injury to the patient or death.

If a user and/or a patient comes to know a serious incident that has occurred in relation to the guide wire, report the fact to the manufacturer and the competent authority of the Member State.

These Instructions for Use apply to the ASAHI Conquest Pro and/or Confianza Pro. For details (length of the guide wire, length of radiopaque section, etc.), refer to the product label.

Descriptions

The basic structure of ASAHI Conquest Pro/Confianza Pro series consists of a taper core wire and coil. The taper core wire and the coil are soldered. The length of the coil is radiopaque so that its position can be easily confirmed under radioscropy.

The distal surface is coated with hydrophilic coating, and the proximal surface is coated with PTFE (polytetrafluoroethylene). At the proximal end of the devices other than those with total lengths of 300cm, there is a connection for an extension wire, which is used to replace the combined device, and ASAHI INTECC's Extension Wire can be connected.

Refer to the table below for the tip load* of this product.

Table 1. Product Specifications

Product Name	Outer Diameter(Distal/Proximal)	Tip Load*
ASAHI Conquest Pro/ Confianza Pro	0.23mm/0.36mm(0.009inch/0.014 inch)	9.0gf
ASAHI Conquest Pro 12/ Confianza Pro 12	0.23mm/0.36mm(0.009inch/0.014 inch)	12.0gf

*Tip load values represent the mean values of product specification values measured according to the internal standard testing method assuming the force applied to the lesion by the guide wire tip.

Indications for Use

This guide wire is intended for general percutaneous coronary intervention (PCI). This guide wire is intended to facilitate the placement of balloon dilatation catheters or other interventional therapeutic devices within cardiovascularure.

This guide wire is not intended for use in the neurovasculature and peripheral vasculature.

Target Treated Population

ASAHI Conquest Pro/Confianza Pro is intended for patients who meet indications for PCI procedures. It is not indicated for use on neonates, infants, or children as well as patients who are pregnant, have impaired renal function, and who have an allergic reaction to contrast media or polyethylene glycol (PEG).

Contraindications

- Confirm the compatibility of guide wire diameter with other interventional device before actual use.
- Do not use the guide wire in patients who have an allergic reaction to contrast media or who have impaired renal function.
- This guide wire contains PEG400 in its hydrophilic coating. Do not use the guide wire in patients who have allergic reactions to polyethylene glycol (PEG).
- Do not use the guide wire for pregnant patients or patients who may be pregnant (X-ray may affect the fetus).
- This guide wire must be used in an institution where emergency surgical operation can be performed immediately. If an emergency surgical operation is unavailable, in the worst case, life-threatening events may occur.
- This guide wire is presterilized with ethylene oxide gas (EOG) and is intended for single use only. Do not reuse or resterilize. If reused or resterilized, the performance or quality of this guide wire and package may be compromised and there is a risk of complications, including infection.

Warnings

- This guide wire must be used only by a physician who is fully trained in PCI treatment.
- Never use this guide wire in combination with catheters having metal parts that may directly contact the surface of this guide wire (e.g. atherectomy catheters, metal dilators), metal needles or metal sheath inserters.
- During use of the guide wire, check and monitor the movement of its tip under fluoroscopy. Do not use in areas of vessel that are not or cannot be visualized.
- Do not connect this guide wire with detachable extension wires produced by manufacturers other than ASAHI INTECC. Do not manipulate this guide wire while the detachable extension wire is connected. Otherwise, the guide wire may be damaged, or the detachable extension wire may be unintentionally detached.
- Do not modify this guide wire for any reason.

Precautions

- Do not use the guide wire after the expiration date indicated on the label. Discard any guide wire that exceeds the expiration date.
- Guide wires are delicate instruments and should be handled carefully.
- If the package is opened or damaged, do not use the guide wire. Do not open the package until just prior to use. Use aseptic technique in handling and using the guide wire.
- Always keep the guide wire in a wet state to prevent damage to its hydrophilic coating.
- Do not wipe this guide wire using an organic solution such as alcohol.
- Take preventive measures against infection after use. Discard this guide wire and package as medical waste.
- Never push, auger, withdraw, or torque this guide wire that meets resistance. Torquing or pushing this guide wire against resistance may cause damage and/or tip separation of this guide wire or direct damage to a vessel. Failure to do so may cause damage to this guide wire. Determine the cause of resistance under fluoroscopy and take any necessary remedial action.
- Always advance or withdraw the guidewires slowly and carefully.
- When torquing this guide wire inside the blood vessel, do not torque continuously in the same direction. When torquing the guide wire, rotate it clockwise and counterclockwise alternately. Do not rotate more than two rotations (720°) in the same direction.
- Flush the guide wire with heparinized and sterilized saline or other suitable solution while removing and reinserting it to prevent air from entering the interventional device. Perform exchange of this guide wire carefully to prevent air entry and/or trauma. When reintroducing the guide wire, confirm that the interventional device tip is free within the vessel lumen and is not against the vessel wall. Failure to do so may result in trauma. Use the radiopaque marker of the interventional device to confirm position.

- Free movement of the guide wire within the interventional device is an important feature of a steerable guide wire system because it gives the user valuable tactile information. Test the system for any resistance prior to use. Adjust or replace the hemostasis valve with an adjustable valve if it is found to inhibit the guide wire movement.
- Do not perform stent placement using more than one guide wire or wire operation through stent strut. Otherwise, the stent may be damaged or the guide wire may break or break apart.
- Before use, confirm that this guide wire is appropriate in size and shape for the procedure and is compatible with the device to be used in combination.
- Refer to the Instruction for Use attached to the medical device (e.g. balloon dilatation catheters or other interventional therapeutic devices) to be used in combination with this guide wire.
- When connecting the detachable extension wire with this guide wire, it should be securely inserted to the boundary line between non PTFE coating and PTFE coating at the proximal end of the guide wire. If any abnormality is found in the connection, stop using the guide wire immediately.
- Inspect the guide wire carefully for bends, kinks, or other damage prior to use and whenever possible during the procedure. Do not use the guide wire if it is damaged.
- When shaping the tip, be careful not to damage the coil.
- Before inserting the guide wire, check which end is the distal end and make sure to insert it from the distal end.
- Maintain a continuous saline flush between the guiding catheter and the interventional device and between the interventional device and the guide wire during the procedure. Flushing prevents contrast media formation and/or clotting on the guide wire and in the catheter lumen.

- entire device. Take note that heparinized and sterilized saline may spill out of the holder tube during at this time.
- c) Fill the lumen of the interventional device with the heparinized and sterilized saline before

- inserting the guide wire.
- d) Release the guide wire from the tail clip and slide the guide wire slowly, push distal end of the guide wire through the holder tube(B).(Figure 1)

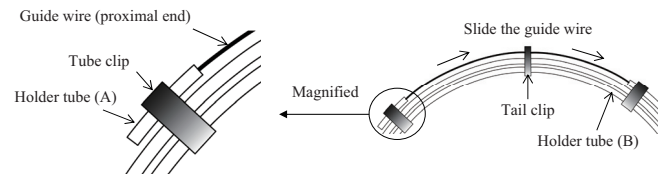
Figure 1

- e) When the distal end of the guide wire is extended 5 to 6 cm beyond the other end of the holder tube, if necessary, shape the tip in accordance with standard practice.
- f) Gently grasp the guide wire which came out from the distal end of the holder tube, at the point as close to the holder tube as possible and pull the guide wire out slowly and carefully.
- g) If resistance is felt when removing the guide wire from the holder tube, inject heparinized and sterilized saline into the holder tube again. Continue removing the guide wire after the resistance is not felt anymore. Be sure to inject enough heparinized and sterilized saline to the holder tube because hydrophilic coating may be damaged when removing the coated guide wire forcefully.
- h) After pulling the guide wire out of the holder tube, inspect it to make sure that it is not damaged.

the guide wire system is just proximal to the tip of the guiding catheter.

- e) Tighten the hemostasis valve of Yd connector to create a seal around the interventional device. Ensure the guide wire movement is still permitted.
- f) Check to make sure the guide wire moves smoothly.
- g) Attach the torque device to the guide wire if necessary.
- h) Advance this guide wire under fluoroscopy to pass through the lesion. Confirm by angiography that the guide wire has passed the narrowed target lesion.
- i) Observe the guide wire movement in the vessels.
- j) Advance the interventional device until the lesion is reached while preventing the guide wire from moving. Ensure that guide wire distal tip and its location in the vessel are visible during interventional device manipulations.

3. Procedures for insertion

- Over-the-wire system
 - Insert the distal end of the guide wire carefully into the guide wire lumen of the interventional device.
 - Advance this guide wire carefully until its tip is just proximal to the interventional device tip.
 - Engage the guiding catheter and insert the interventional device system (with guide wire) into the Y connector.
 - Advance the interventional device system through the guiding catheter until the tip of

■ Rapid exchange system

- Engage the guiding catheter.
- Insert the guide wire introducer into the Y connector of the guiding catheter.
- Carefully insert the guide wire tip into the introducer.
- Advance the guide wire through the guiding catheter under fluoroscopy until the guide wire tip is just proximal to the tip of the guiding catheter.
- Attach the torque device to the guide wire if necessary.

Malfunctions and Adverse Events

Possible malfunctions and adverse effects through usage of this guide wire include, but are not limited to the following:

1. Malfunctions

Failure to cross a lesion	Breakage or bending of the guide wire	Damage such as separation
Difficulty in removal	Debris from a damaged guide wire	

2. Adverse events

Death	Infection	Vessel dissection
Vessel perforation	Bleeding complications	Embolism
Thrombus	Infarction	Ischemia
Arrhythmia	Vasospasm	Vascular occlusion
Blood pressure reduction	Allergic reaction	

How to Use

1. Inspection prior to use

- Before use, inspect carefully and confirm the guide wire and package are undamaged.
- Before use, confirm that the guide wire is compatible with the interventional devices to be used.

2. Preparation

- Select the most suitable guide wire for the affected area and remove the holder tube containing the guide wire from the sterile pack.
- Fill up the holder tube with heparinized and sterilized saline using a syringe to soak the

- f) Advance the guide wire under fluoroscopy to pass through the lesion. Confirm by angiography that the guide wire has passed the narrowed target lesion.
 - g) Observe the guide wire movement in the vessels.
 - h) Remove the guide wire torque device and the guide wire introducer.
 - i) Track the interventional device over the guide wire while preventing the guide wire from moving, and advance until the lesion is reached. Ensure that the guide wire distal tip and its location in the vessel are visible during interventional device manipulations.
4. Procedures to exchange the guide wire
- Over-the-wire system
 - a) Remove this guide wire slowly while monitoring the movement of this guide wire under fluoroscopy.
 - b) Insert the next guide wire in accordance with the directions in this “How to Use” section.

Storage Method

Do not keep the guide wire and package in a bent and/or heavily-loaded condition. This guide wire and package must be kept out of water. Store in a cool, dark, and dry place.

Expiration Date

The expiration date is indicated on the label of the guide wire package.

Contents

5 pieces / box

Liability Disclaimer

By no means shall “ASAHI INTECC CO., LTD. and its affiliated companies” (hereinafter referred to as the “Company”) be liable for accidents, personal injuries, adverse events due to any improper use of the product(s) or any other use inconsistent with these instructions. In no event shall the Company be liable for any damages either (i) arising out of storage of the product(s) after the shipment from the Company or (ii) due to selection of patients, surgery techniques, or any other medical activities by the medical institution that uses the product(s).

“ASAHI Conquest” and “Confianza Pro” are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

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©2024 ASAHI INTECC CO., LTD.
Emergo Europe B.V.