


Refinity is used with the SpinVision automatic pullback device

- 15 cm, automatic pullbacks for accurate vessel and lesion measurements
- Manual pullback is available when desired

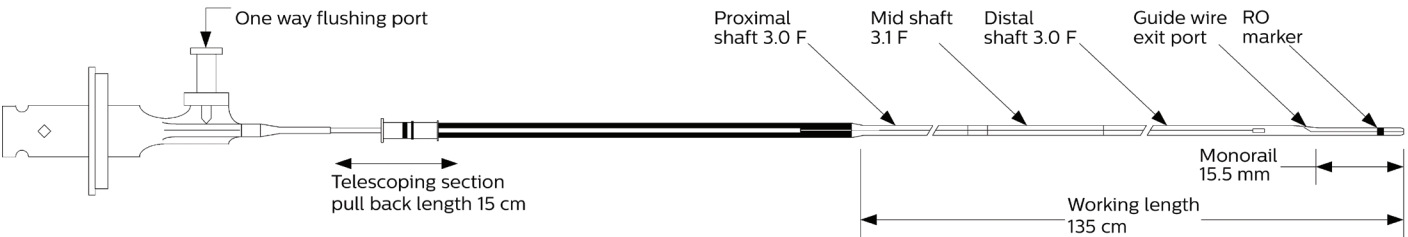

Dimensional attribute	Specification
Imaging core pullback length	15 cm
Catheter working length	135 cm
Proximal shaft profile	3.0 F (0.39")
Guide wire compatibility	0.014"
Guide catheter compatibility	5 F (ID ≥ .056")
Imaging window profile	3.0 F (.39")
Monorail length	15.5 mm
Transducer to tip length	20.5 mm
RO marker to tip length	5 mm

Ordering information	
Refinity short-tip catheter	89900
SpinVision option	891 SV
SpinVision sterile drapes (order by each)	891 BG

Refinity ST IVUS Catheter

Discover more choices for your complex PCI procedures

PHILIPS



With a low crossing profile and high definition imaging, Refinity ST is designed to meet the needs of your most complex PCI cases.


Features designed for enhanced deliverability

- Large 0.026 inch Asahi Intecc drive cable for enhanced pushability
- Smooth proximal to midshaft transition
- GlyDx hydrophilic coating for enhanced pushability and trackability
- Soft and flexible tip facilitates distal access



45 MHz high definition image

- Clear images for accurate planning, delivery and optimization of coronary stents

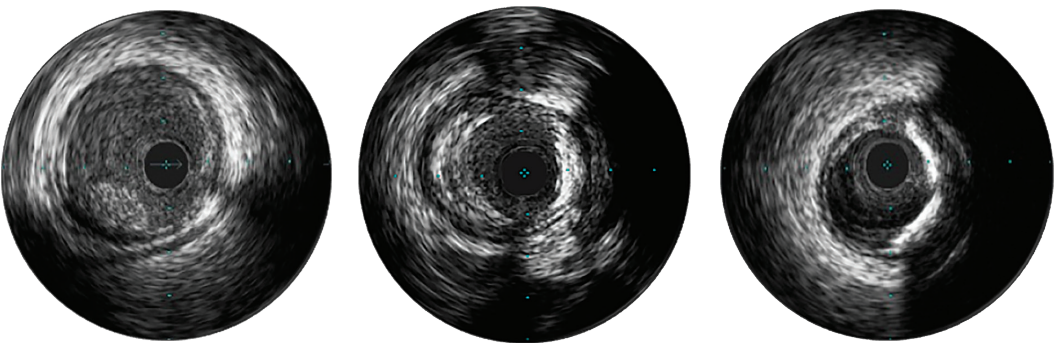

Similar catheters

	Refinity	OptiCross HD
Frequency/type	45 MHz Rotational IVUS	60 MHz Rotational IVUS
Crossing profile	3.0 F	3.1 F
Imaging window profile	3.0 F	2.6 F
Proximal shaft profile	3.0 F	3.0 F
Tip entry profile	0.025"	0.026"
Minimum guide catheter	5F	5F
Pullback length	150 mm	100 mm

More options for your patients

Physician preferences and challenging lesions require options. **Only Philips** provides a wide range of solutions from a single system.

- Digital IVUS
- Rotational IVUS
- Peripheral IVUS
- iFR and FFR modalities

